

Single Crystal Analysis by XtaLAB Synergy-ED

Product used : Electron Diffractometer


A novel method for determining the structure of submicron particles



Left : Backscattered electron image of Diflubenzuron in LV-mode, JCM-7000 NeoScope™


Right : Electron diffraction structure analysis of Diflubenzuron, XtaLAB Synergy-ED

XtaLAB Synergy-ED is a fully optimized electron diffractometer for microcrystals analysis, jointly developed by both companies: Rigaku Corporation and JEOL Ltd. The XtaLAB Synergy-ED features Rigaku's high-speed and high-sensitivity photon-counting detector (HyPix-ED) and a state-of-the-art instrument control with a single crystal analysis software platform (CrysAlis^{Pro} for ED). The key characteristic of this product provides researchers an easy and efficient platform for electron crystallography. In the example above, the right-hand shows the result of structural analysis of Diflubenzuron (C₁₄H₉ClF₂N₂O₂).

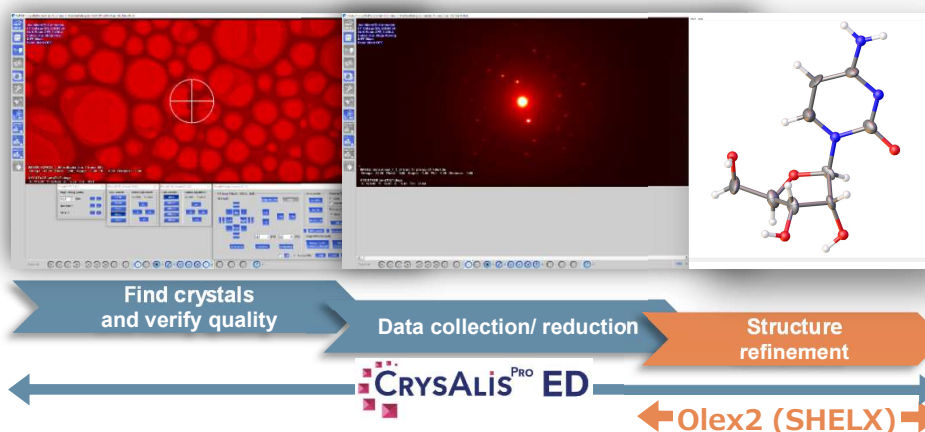


Sample preparation

No sample preparation is required: load your sample directly onto a TEM grid as follows. If your particles are large, grind them with a microscope slide or spatula to a fine particle size.

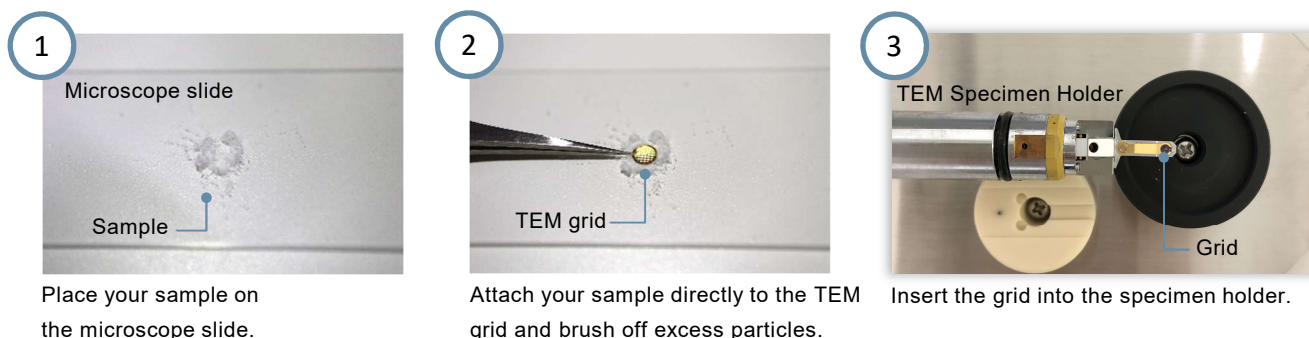

3D ED/microED experiments

CrysAlis^{Pro} for ED allows easy switching between the visual mode and the diffraction mode. In the visual mode, click on an appropriate particle and the Diff.exp button in the MicroED Stage Control menu to switch into the diffraction mode. And then click the Start experiment button in the MicroED/3D ED DC window to start the electron diffraction observation automatically. The rotation range used for a typical measurement is -60 to +60 degrees.

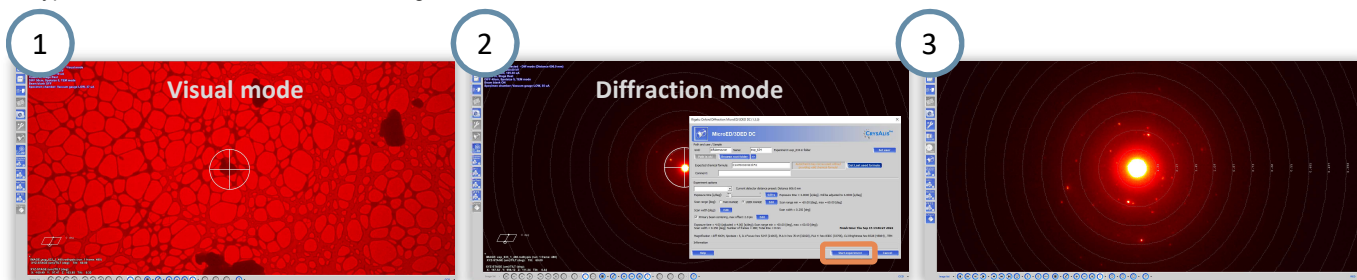


1 Search for crystals in the visual mode.

2 Click the Start experiment button.

3 Start the electron diffraction observation automatically.

Structure analysis

In CrysAlis^{Pro} for ED, Autochem performs fully automatic structure determination and refinement during data collection. In addition, it is also possible to manually analyze the structure with efficient and effective functions and merge the dataset obtained from multiple crystals to improve the completeness of the data.

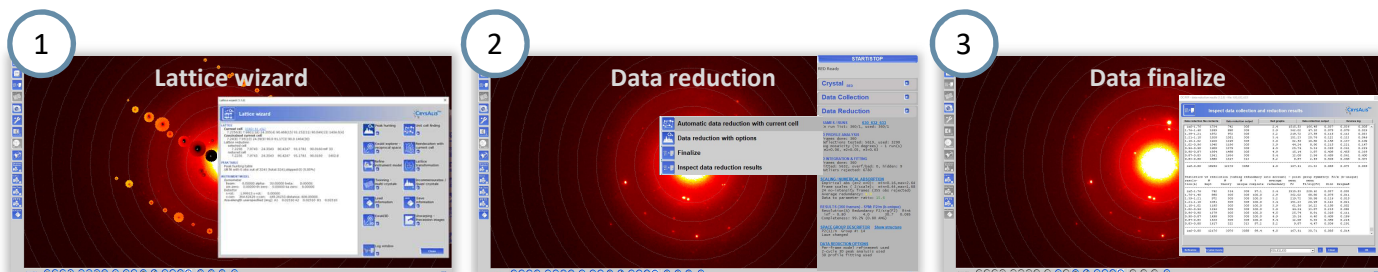


1 Peak detection and unit cell determination.

2 Automatic or interactive data reduction.

3 Data finalize

Structure refinement

CrysAlis^{Pro} for ED allows seamlessly switching to Olex2. The intuitive user interface allows you to build and refine your molecular structure.

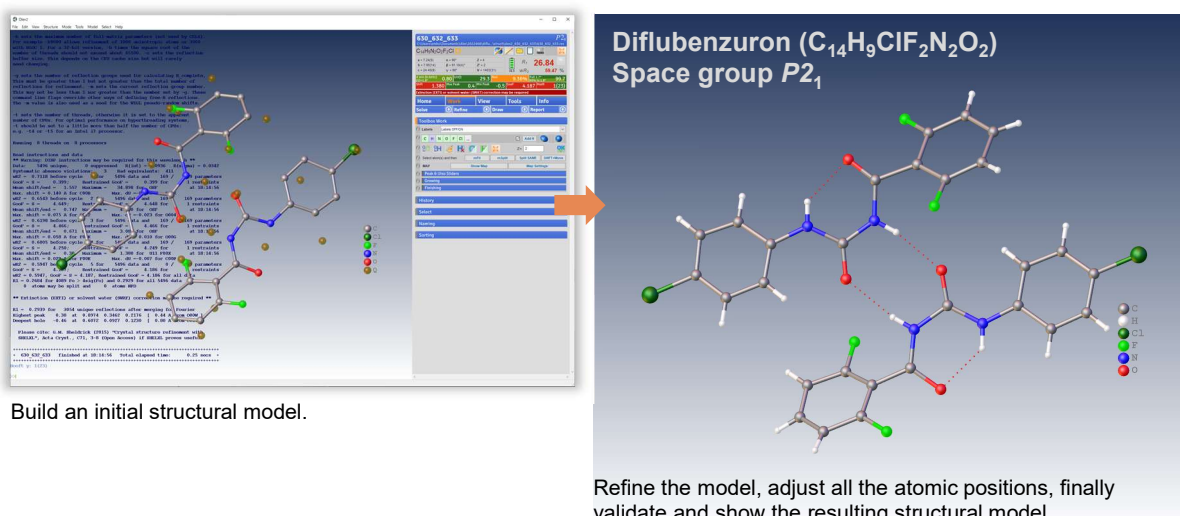

Build an initial structural model.

Refine the model, adjust all the atomic positions, finally validate and show the resulting structural model.

Copyright © 2022 JEOL Ltd.

Certain products in this brochure are controlled under the "Foreign Exchange and Foreign Trade Law" of Japan in compliance with international security export control. JEOL Ltd. must provide the Japanese Government with "End-user's Statement of Assurance" and "End-use Certificate" in order to obtain the export license needed for export from Japan. If the product to be exported is in this category, the end user will be asked to fill in these certificate forms.



3-1-2 Musashino Akishima Tokyo 196-8558 Japan Sales Division Tel. +81-3-6262-3560 Fax. +81-3-6262-3577
www.jeol.com ISO 9001 • ISO 14001 Certified

• AUSTRALIA & NEW ZEALAND • BELGIUM • BRAZIL • CANADA • CHINA • EGYPT • FRANCE • GERMANY • GREAT BRITAIN & IRELAND • ITALY • KOREA • MALAYSIA • MEXICO • RUSSIA • SCANDINAVIA
• SINGAPORE • TAIWAN • THE NETHERLANDS • USA

