

## Structural Analysis of Polyethylene Terephthalates with Different Crystallinity using Reactive Pyrolysis GC-TOFMS and NMR

Product used : Mass Spectrometer (MS), Nuclear Magnetic Resonance (NMR)

Polyethylene terephthalate (PET) is a thermoplastic polyester that is obtained by polycondensation of ethylene glycol and terephthalic acid and has excellent transparency, toughness, rigidity, and heat resistance (Figure 1). PET can be roughly divided into two types, crystalline PET (C-PET) and amorphous PET (A-PET), depending on the processing method. C-PET has the characteristics of high strength and heat resistance due to the regular arrangement of the molecules in the crystallized part, which increases the density. A-PET is characterized by high impact strength and easy bending. However, the amorphous portion of A-PET slowly crystallizes due to long-term use and heat history, causing changes over time. Changes in density generate internal stress that cuts polymer chains, which can lead to deterioration in flexibility, impact resistance, strength, etc. Therefore, a polymer was devised in which about 30-40% of the ethylene glycol in PET was replaced with cyclohexanedimethanol, and this is called glycol-modified PET (G-PET, PETG) (Figure 2). G-PET is treated as an amorphous resin because the polymer does not crystallize even during the molding process.

In this application note, two types of commercially available PET resins are analyzed by reactive pyrolysis GC-TOFMS and NMR, and the results of analysis to determine whether they are PET or G-PET are reported.

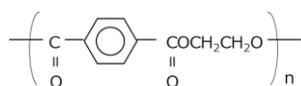


Figure 1 PET (A-PET & C-PET) structural formula

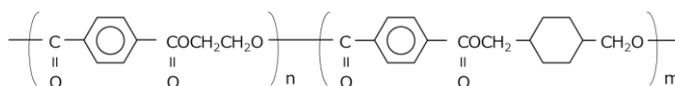


Figure 2 G-PET structural formula

### Experiments

Both commercially available PET film and PET plate were freeze-ground and used as samples. The measurement conditions for each analytical instrument are shown below.

#### Reactive pyrolysis GC-TOFMS

A gas chromatograph time-of-flight mass spectrometer JMS-T2000GC (JEOL Ltd.) equipped with a pyrolyzer EGA/PY-3030D (Frontier Labs Ltd.) was used for the measurements. Samples were subjected to measurements together with 10  $\mu$ L of tetramethylammonium hydroxide (TMAH) 25% (w/w) methanol solution (EI method:  $n=3$ , FI method:  $n=1$ ). The sample amount was 0.2 mg for the EI method and 0.5 mg for the FI method. Using msFineAnalysis (JEOL Ltd.), the obtained data were subjected to differential analysis between the two samples of PET film and PET plate, followed by qualitative analysis and structural elucidation of the components characteristic of G-PET. Other detailed conditions are shown in Table 1.

#### NMR

The sample was dissolved in deuterated trifluoroacetic acid (TFA-d1) and  $^1\text{H}$  NMR measurements were performed using JNM-ECZL400S (JEOL Ltd.). The NMR database of PolyInfo (NIMS) [1] was used as a reference for the structural analysis.



Table 1 Py-GC-TOFMS measurement and analysis conditions

Pyrolysis conditions		MS conditions	
Pyrolyzer	EGA/PY-3030D (Frontier Labs Ltd.)	Spectrometer	JMS-T2000GC (JEOL Ltd.)
Pyrolysis Temperature	400°C	Ion Source	EI/FI combination ion source
GC conditions		Ionization	EI+:70 eV, 300 $\mu$ A
Column	DB-5MS UI (Agilent Technologies)		FI+: -10 kV
	30 m x 0.25 mm, 0.25 $\mu$ m	Mass Range	$m/z$ 29 - 600
Oven Temperature	40°C (2 min) - 20°C/min		
	- 320°C (10 min)		
Injection Mode	Split mode (100:1)		
Carrier flow	He: 1.0 mL/min		

## Results

### Measurement results of reactive pyrolysis GC-TOFMS

Figure 3 shows EI TICC of reactive pyrolysis GC-TOFMS (blue: PET film, red: PET plate). Peaks 1 to 3 were observed as common components for both samples. In addition, peaks 4 to 11 were observed as components characteristic of the PET plate. Table 2 shows the integrated analysis results by msFineAnalysis of peaks 1 to 11 from the PET plate. The integrated analysis is a method to uniquely determine the molecular formula by integrating and analyzing the library database (DB) search results using the EI spectrum and the composition estimation results of the molecular ions in the FI spectrum. From this, it was found that peaks 1 to 3 are methylated derivatives of terephthalic acid and ethylene glycol, which are reaction pyrolysis products of PET. From the results of NIST library database (DB) search and molecular ion composition estimation of the FI spectrum, peaks 6 and 7 were presumed to be structure B or its isomer, and peaks 8 and 9 were presumed to be structure C or its isomer, respectively. Although peaks 4 and 5 could not be found in the DB search, they could be assumed to be structure A or its isomers from the results of EI spectrum and molecular ion composition estimation of FI spectrum. As described above, peaks 4 to 9 were found to be three types of reaction pyrolysis products derived from cyclohexanedimethanol and having different numbers of methyl groups. It was also speculated that there are two isomers for each type.

Next, peak 10 was presumed to be structure D from the DB search results and molecular ion composition estimation results of the FI spectrum, and was a reaction pyrolysis product reflecting the structure of PET. On the other hand, there was no corresponding compound for peak 11 in the DB search. Structural analysis was carried out from characteristic fragment ions ( $m/z$  93, 108, 126) when comparing the EI spectrum with that of peak 10 and the result of molecular ion composition estimation of the FI spectrum. As a result, it was presumed to be structure E in which terephthalic acid and cyclohexanedimethanol were bound.

From the above results, it was presumed that the PET film was PET (A-PET or C-PET) and the PET plate was G-PET.

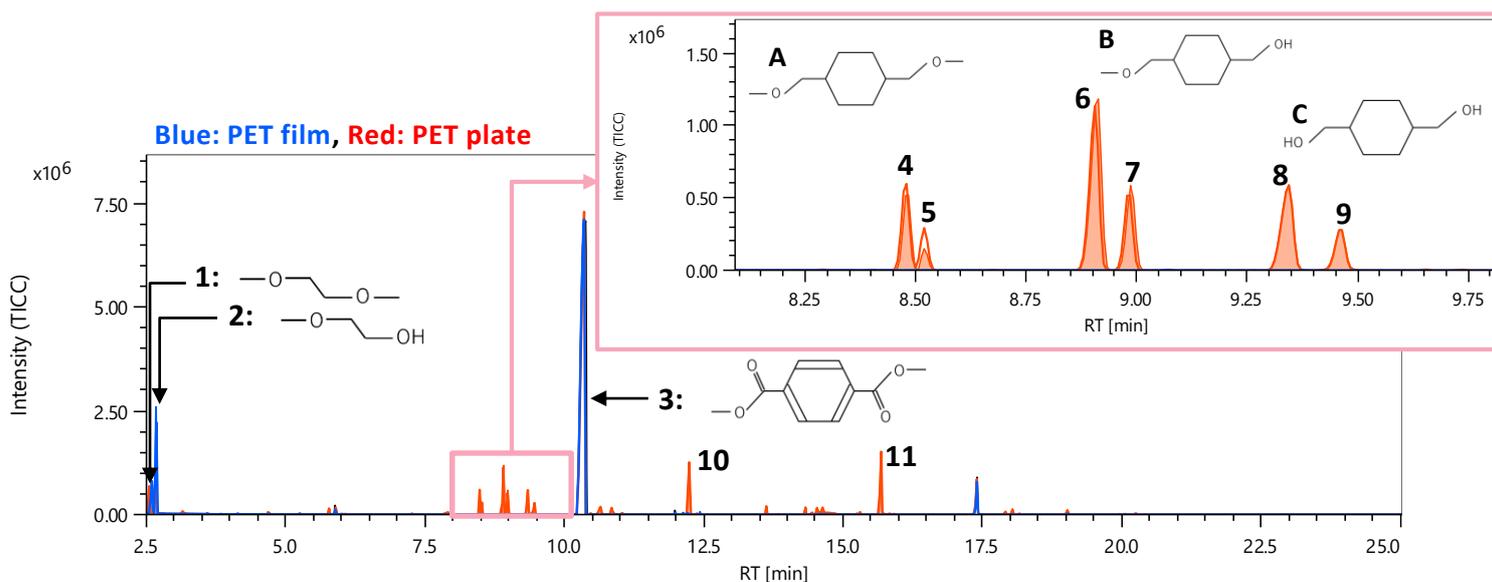


Figure 3 TICC of EI data

Table 2 Integrated qualitative analysis result of PET plate

ID	General			Total Result							
	RT [min]	Height [%]	IM $m/z$	Library Name	Lib.	Similarity	Formula	DBE	Calculated $m/z$	Mass Error [mDa]	EI Fragment Coverage
1	2.55	9.40	76.05187	Ethanol, 2-methoxy-	mainlib	910	C3 H8 O2	0.0	76.05188	-0.01	100
2	2.60	10.35	90.06737	Ethane, 1,2-dimethoxy-	mainlib	771	C4 H10 O2	0.0	90.06753	-0.16	88
3	10.35	100.00	194.05831	1,4-Benzenedicarboxylic acid, dimethyl ester	mainlib	932	C10 H10 O4	6.0	194.05736	0.95	100
4	8.48	7.34	172.14621	-	-	-	C10 H20 O2	1.0	172.14578	0.43	100
5	8.52	2.11	172.14641	-	-	-	C10 H20 O2	1.0	172.14578	0.63	100
6	8.91	16.74	158.13057	cyclohexanemethanol, 4-(methoxymethyl)-	mainlib	912	C9 H18 O2	1.0	158.13013	0.44	100
7	8.98	8.29	158.13052	cyclohexanemethanol, 4-(methoxymethyl)-	mainlib	844	C9 H18 O2	1.0	158.13013	0.39	100
8	9.34	8.32	144.11460	1,4-Cyclohexanedimethanol, trans-	mainlib	973	C8 H16 O2	1.0	144.11448	0.12	100
9	9.46	3.94	144.11466	1,4-Cyclohexanedimethanol, trans-	mainlib	856	C8 H16 O2	1.0	144.11448	0.18	100
10	12.44	1.02	266.07821	1,2-Benzenedicarboxylic acid, 2-ethoxy-2-oxoethyl methyl ester	mainlib	712	C13 H14 O6	7.0	266.07849	-0.28	100
11	15.69	21.42	306.14702	-	-	-	C17 H22 O5	7.0	306.14618	0.84	100

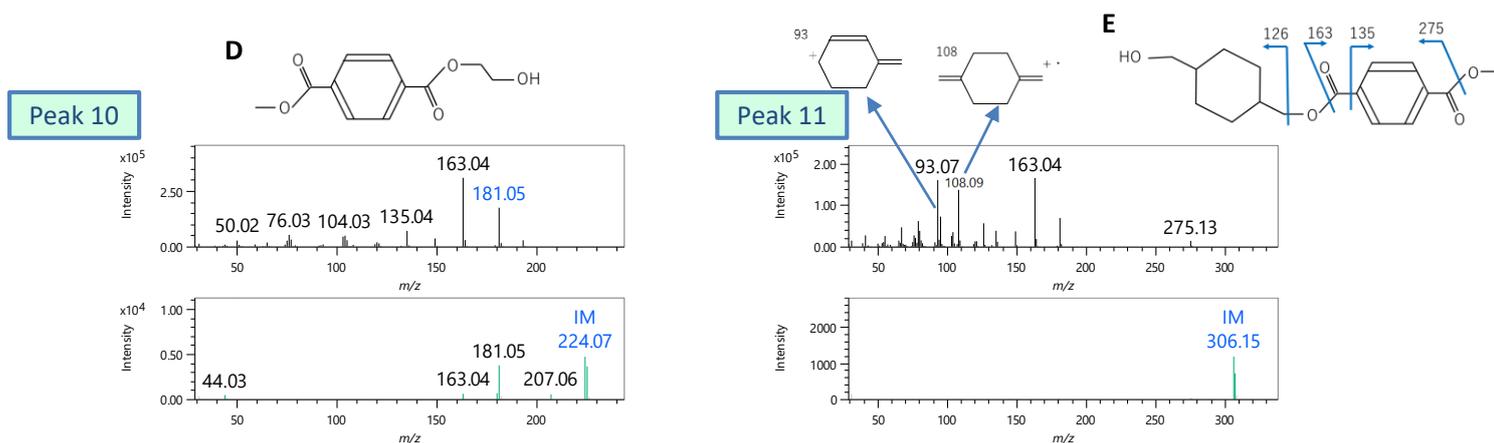


Figure 4 Mass spectra of peaks 10 and 11

## <sup>1</sup>H NMR measurement results

Figure 5 shows the <sup>1</sup>H NMR spectra of the PET film and the PET plate. Peaks 1 and 2 derived from terephthalic acid and ethylene glycol were commonly observed in the spectra of both samples. On the other hand, multiple peaks (peaks 3 to 12) were observed around chemical shifts of 1 to 5 ppm only for the PET plate. The PET plate was presumed to be G-PET from the results of reactive pyrolysis GC-TOFMS, so we searched PoLyInfo by polymer name and compared the NMR spectra listed with the observed NMR spectrum, and found high similarity. After confirming that, peak assignment was performed. As a result, peaks 3 to 12 were found to be derived from cyclohexanedimethanol. In the NMR spectrum, the peaks derived from the structural isomers of cyclohexanedimethanol are separated, and it is characteristic that the abundance ratio of the isomers can be discussed. Also, when the ratio of ethylene glycol and cyclohexanedimethanol was calculated from each peak area value for the PET plate, the ratio was approximately 2:1.

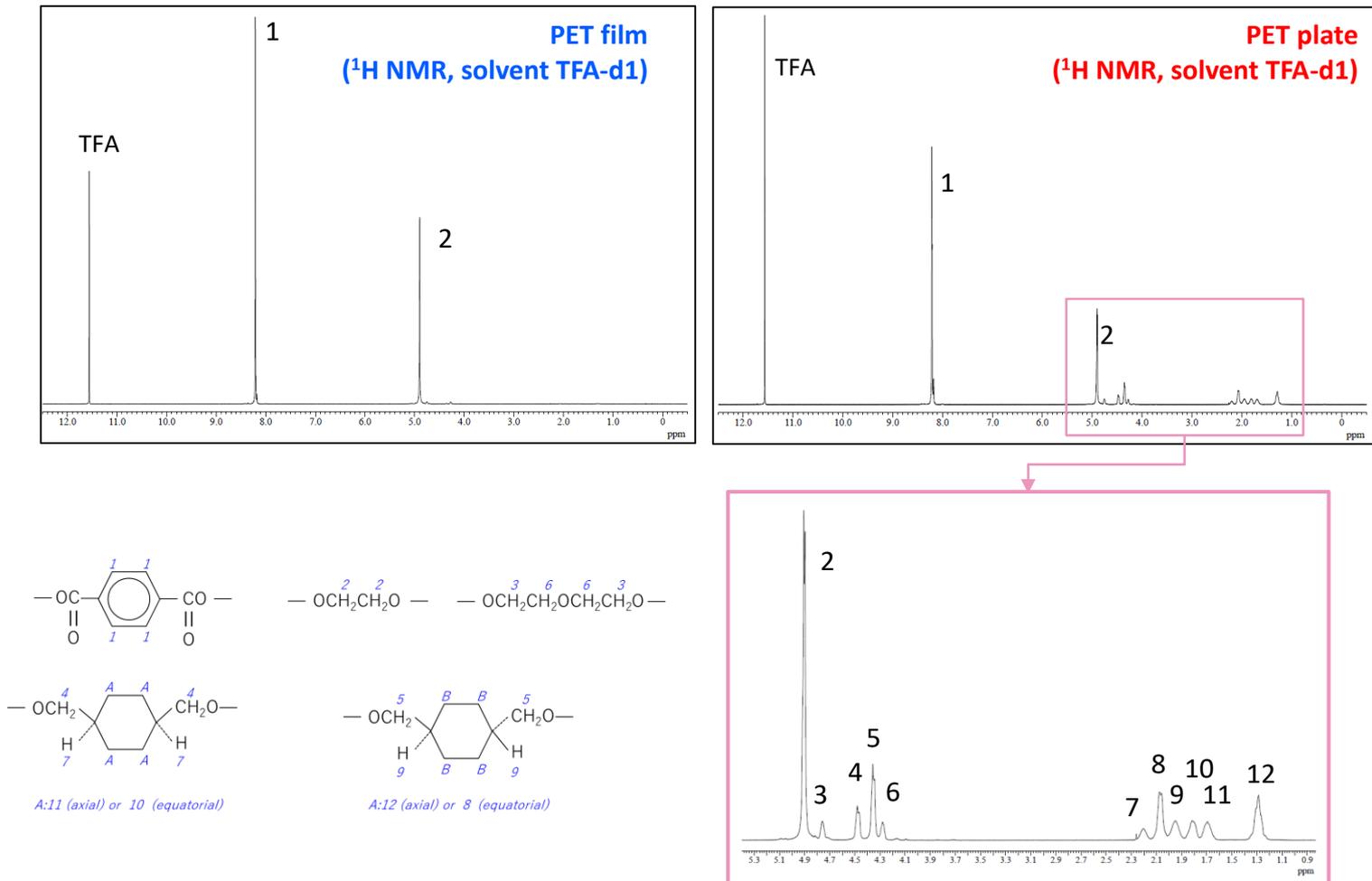


Figure 5 <sup>1</sup>H NMR spectra of PET film and PET plate

## Summary

In this application note, we reported the results of structural analysis of commercially available PET film and PET plate by reactive pyrolysis GC-TOFMS and <sup>1</sup>H NMR. Reactive pyrolysis GC-TOFMS clarified the difference in structure between the PET film and the PET plate, and was able to deduce that they were PET (A-PET or C-PET) and G-PET, respectively. Based on this information, the analysis of the <sup>1</sup>H NMR results was facilitated, and it was found that the structural isomers and the ratio of ethylene glycol/cyclohexanedimethanol could be discussed. In this way, it is effective to comprehensively utilize the information obtained from each technique in structural analysis of polymers.

## Reference

- [1] NMR Database, <https://polymer.nims.go.jp> [accessed 2022-08-31]

Copyright © 2023 JEOL Ltd.

Certain products in this brochure are controlled under the "Foreign Exchange and Foreign Trade Law" of Japan in compliance with international security export control. JEOL Ltd. must provide the Japanese Government with "End-user's Statement of Assurance" and "End-use Certificate" in order to obtain the export license needed for export from Japan. If the product to be exported is in this category, the end user will be asked to fill in these certificate forms.

