


JEOL Products Guide

Products Guide 2026

Nanometrology

Fabrication

Characterization

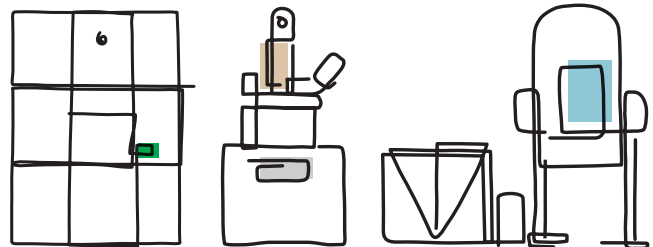

Introduction

Thank you for your interest in JEOL products and services.

JEOL designs and manufactures scientific instruments for high-level research and development activities. Our customers include scientists and engineers working in leading-edge academic and industrial laboratories around the world. JEOL products and services enable them to pursue a variety of R&D applications that require high resolution imaging and analytical capabilities such as: basic observation and analysis, environmental science, information technology, semiconductor production, biotechnology, nanotechnology, and a broad range of industrial endeavors.

Experts involved in the studies of medicine, biology, biochemistry, agriculture, materials science, metallurgy, ceramics, chemistry, petroleum, pharmacy, semiconductors and electronic materials have been using JEOL products for more than 75 years. Our new products are easier-to-use than ever before and contribute to a high level of quality assurance and quality control during the production process.

This Product Guide presents the most current high-performance solutions from JEOL to meet your R&D requirements for electron optics, analytical, semiconductor, and industrial instruments and equipment. For more details or information about any of our products, please contact your nearest JEOL office.

JEOL Ltd.

Company Profile

Name: JEOL Ltd.

Address (head office): 1-2, Musashino 3-chome, Akishima, Tokyo 196-8558,


Japan Establishment: May 30, 1949

Capital: 21,394 million yen (as of the end of March 2025)

Consolidated sales: 196,695 million yen (as of the end of March 2025)

Number of employees: 3,604 (for JEOL Group as of the end of March 2025)

Contents

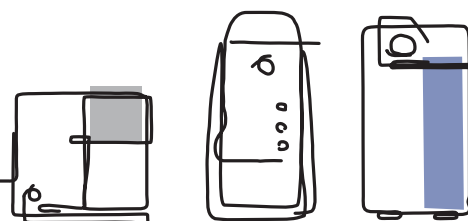

Scientific/Metrology Instruments	■ Transmission Electron Microscopes (TEM)	P. 3
	■ Scanning Electron Microscopes (SEM)	P. 5
	■ Ion Beam Application Equipment	P. 7
	■ Peripheral Equipment	P. 8
	■ Instruments for Microarea and Surface Analysis	P. 9
	■ X-ray Fluorescence Spectrometers	P. 10
	■ Electron Diffractometer	P. 10
	■ Mass Spectrometers	P. 11
	■ NMR Spectrometers	P. 13
	■ ESR Spectrometers	P. 14
Industrial Equipment	■ Semiconductor Equipment	P. 15
	■ Metal AM Machine	P. 16
	■ Industrial Equipment for thin-film formation and material processing	P. 16

* Some instrument photographs include optional attachments.

* Specifications subject to change without notice.

* This catalog includes products not offered in some territories.

* Windows is a registered trademark of Microsoft Corp. in the U.S. and other countries; Macintosh is trademark of Apple Computer Inc., registered in the U.S. and other countries; UNIX is a registered trademark of X/Open Corp. of the USA; NetWare is a registered trademark of Novell Inc. in the U.S. and other countries; NetWare Lite is a trademark of Novell Inc. in the U.S. and other countries; Pentium is a registered trademark of Intel Corp. of the USA.



Transmission Electron Microscopes



GRAND ARM™ 3

Atomic Resolution Analytical Electron Microscope

JEOL's flagship model, 300 kV atomic resolution GRAND ARM™ 3, comes newly equipped with enhanced observation and analysis capabilities. JEOL's joint development with CEOS GmbH, Germany has enabled the GRAND ARM™ 3 to be outfitted with the latest, most advanced aberration correctors, Large Aperture STEM Corrector (LASCOR) and Advanced TEM Corrector (ATCOR). These correctors allow for imaging and analysis under much larger convergence angles than previously attainable; as a result, the GRAND ARM™ 3 yields unprecedented ultra-high resolution TEM and STEM results.

- STEM Resolution: 0.049 nm (FHP2, 300 kV, with STEM Corrector)
- TEM Resolution: 0.050 nm (FHP2, 300 kV, with TEM Corrector)
- Accelerating voltage: 40 to 300 kV (Standard 80 kV, 300 kV)
- Electron gun: Cold field emission gun (CFEG)



NEOARM

Atomic Resolution Analytical Electron Microscope

The JEM-ARM200F "NEOARM" comes with JEOL's unique cold field emission gun (Cold-FEG) and a new Cs corrector (ASCOR) that compensates for higher order aberrations. The combination of a Cold-FEG and ASCOR enables atomic-resolution imaging at not only 200 kV accelerating voltage, but also a low voltage of 30 kV.

"NEOARM" is also equipped with an automated aberration correction system. This system incorporates JEOL's new aberration correction algorithm for automatic fast and precise aberration correction. JEOL delivers high-throughput atomic resolution images to everyone.

- STEM Resolution: 0.071 nm (UHR, 200 kV, with STEM Corrector)
- TEM Resolution: 0.07 nm (UHR, 200 kV, with TEM Corrector)
- Accelerating voltage: 30 to 200 kV (Standard 80 kV, 200 kV)
- Electron gun: Cold field emission gun (CFEG)



CRYO ARM™ 300S

Field Emission Cryo-Electron Microscope

The "CRYO ARM™ 300S" is a cryo-electron microscope that specializes in the observation of electron beam-sensitive specimens, such as proteins, for single particle analysis, tomography and 3DED/MicroED. This system offers improved stability, throughput and ease of use compared to the previous generation of cryo-EMs. Moreover, this is an all-in-one system that can handle everything from screening to data acquisition, allowing for more flexibility in operation at customer sites to meet the needs of the facility. These improvements allow users to obtain high quality images by simple operation even for those who have never used an electron microscope before.

- Accelerating voltage: 300 kV
- Energy filter: In-column Omega energy filter
- Maximum tilt angle: $\pm 70^\circ$
- Specimen cooling temperature: 105 K or less
- Specimen storage: Up to 12 specimens can be held.
- Electron gun: Cold field emission gun (CFEG)



CRYO ARM™ 200S

Field Emission Cryo-Electron Microscope

The "CRYO ARM™ 200S" is a cryo-electron microscope that specializes in single particle analysis used for 3D structure analysis of proteins, etc. In line with the JEOL CRYO ARM™ series, this system is equipped with a CFEG, an in-column energy filter, a side-entry liquid-nitrogen cooling stage and an automated specimen exchange system. The automated exchange system enables the storing of up to 12 samples and allows for the exchange of one or more samples, thus enabling flexible observation scheduling by multiple users. Moreover, this cryo-EM incorporates a newly-designed high-resolution pole-piece for acquiring higher-quality images.

- Accelerating voltage: 200 kV
- Energy filter: In-column Omega energy filter
- Specimen cooling temperature: 105 K or less
- Specimen storage: Up to 12 specimens can be held.
- Electron gun: Cold field emission gun (CFEG)



ACE eye™

High Throughput Analytical Electron Microscope

The ACE eye™ is a semiconductor-dedicated TEM that delivers fully automated, end-to-end operation—from TEM/STEM imaging to EDS analysis—via a recipe-based workflow. By standardizing measurement conditions, it ensures outstanding reproducibility while offering flexible operation modes, ranging from semi-automated to fully-automated workflows. Equipped with AI-driven automatic aberration correction (ASCOR, as a standard) and advanced automation algorithms, ACE eye™ continuously maintains optimal measurement conditions, enabling highly precise, operator-independent results. In addition, by leveraging proven hardware technologies and a design optimized for long-term stable operation, ACE eye™ delivers the usability and reliability required for semiconductor analysis.

- STEM Resolution: 0.10 nm (200 kV, with STEM Corrector)
- TEM Resolution: 0.10 nm (200 kV)
- Accelerating voltage: 60 to 200 kV (standard 80 kV, 200 kV)
- Electron gun: Cold field emission gun (CFEG)



JEM-F200

Multi-purpose Electron Microscope

The JEM-F200 (F2) is a next-generation, multi-purpose electron microscope designed for today's diverse analytical needs. Featuring a new probe-forming lens system, cold field emission gun (CFEG), and Dual SDD system, it delivers high-resolution imaging and efficient, high-throughput materials analysis. The integrated SPECPORTER™ system enables automated holder insertion and retraction, ensuring smooth and reliable operation.

- STEM Resolution: 0.14 nm (UHR, 200 kV, CFEG)
- TEM Resolution: 0.19 nm (Point image, UHR, 200 kV, CFEG)
0.10 nm (Lattice image, UHR, 200 kV, CFEG)
- Accelerating voltage: 20 to 200 kV (Standard 80 kV, 200 kV)
- Electron gun: Cold field emission gun (CFEG) / Thermal field-emission electron gun (TFEG)



JEM-2100Plus

Electron Microscope

The JEM-2100Plus is a multi-purpose transmission electron microscope, which combines the proven JEM-2100 optic system with an advanced control system for enhanced ease of operation. Achieving superior performance through intuitive operation, the JEM-2100Plus provides solutions to a wide range of applications from materials science to medical/biological studies.

- STEM Resolution: 1.0 nm (UHR, 200 kV)
- TEM Resolution: 0.14 nm (Lattice image, UHR, 200 kV)
- Accelerating voltage: 80 to 200 kV



JEM-120i

Electron Microscope

The JEM-120i is a 120 kV transmission electron microscope designed around the core principles of "Compact", "Easy to Use", and "Expandable". With its modern design and intuitive operation, the JEM-120i streamlines workflows from routine imaging to maintenance—making high-quality TEM analysis accessible to users of all experience levels.

- TEM Resolution: 0.2 nm (Lattice image, HC)/0.14 nm (Lattice image, HR)
- Accelerating voltage: 20 to 120 kV
- Maximum tilt angle (Tilt-X): $\pm 80^\circ$ (HC/HR)

Note: The high tilt specimen holder is required as option.



Scanning Electron Microscopes



JSM-IT810
(SHL version)

JSM-IT810

Schottky Field Emission Scanning Electron Microscope

The JSM-IT810 Series represents the flagship of our SEM product line, enabling acquisition of ultra-high-resolution data with exceptional ease of use. Its superior performance is powered by advanced JEOL technologies, including an in-lens Schottky Plus field emission electron gun, the Neo Engine electron optical control system, and SEM Center instrument control software for seamless operation. The system can be equipped with a fully integrated JEOL energy-dispersive X-ray spectrometer (EDS), providing real-time elemental analysis to complement high-resolution imaging. In addition, no code automation tools are built-in for automating image and EDS analysis providing a streamlined and efficient workflow. Choose between three types of objective lens: a hybrid lens (HL), a super hybrid lens (SHL) and semi-in-lens (SIL). SHL and SIL each have two versions, standard and PRIME.

- Resolution: HL: 0.7 nm (20 kV), 1.3 nm (1 kV)
SHL: 0.6 nm (15 kV), 1.1 nm (1 kV)
SHL PRIME: 0.5 nm (15 kV), 0.7 nm (1 kV)
SIL: 0.6 nm (15 kV), 1.0 nm (1 kV)
SIL PRIME: 0.5 nm (15 kV), 0.7 nm (1 kV)
- Landing voltage: 0.01 to 30 kV
- Direct magnification: HL/SHL: $\times 10$ to $\times 2,000,000$ (128 $\times$ 96 mm print size)
SIL: $\times 25$ to $\times 2,000,000$ (128 $\times$ 96 mm print size)
- Display magnification: HL/SHL: $\times 27$ to $\times 5,480,000$ (1,280 $\times$ 960 pixels on display)
SIL: $\times 69$ to $\times 5,480,000$ (1,280 $\times$ 960 pixels on display)



JSM-IT710HR

Scanning Electron Microscope

- Clear visibility promotes new discovery -

The JSM-IT710HR is the latest HR model, designed for effortless acquisition of high-resolution images. This model features seamless navigation from optical to SEM imaging, a new multi-segmented backscatter electron detector for live 3D imaging, and auto functions from alignment to focus for fast clear and sharp images. A new hybrid low vacuum secondary electron detector expands observation options by providing both a secondary and photon signal. Productivity enhancing no code automation functions are built-in for automation of image collection and includes EDS when JEOL's fully embedded EDS is configured with the system.

- Resolution: High vacuum mode: 1.0 nm (20 kV), 3.0 nm (1 kV)
Low vacuum mode: 4.0 nm (30 kV, BED)
Analysis performance: 3.0 nm (15 kV, probe current 3 nA)
- Landing voltage: 0.5 to 30 kV
- Direct magnification: $\times 5$ to $\times 600,000$
(defined with a print size of 128 $\times$ 96 mm)
- Display magnification: $\times 15$ to $\times 1,767,305$, Image magnification on the monitor
(defined with a display size of 377 $\times$ 283 mm)



JSM-IT510

Scanning Electron Microscope

The JSM-IT510 is the most versatile of our thermionic-emission SEMs and comes with a large multipurpose chamber that can accommodate large specimens (200 mm dia.). This multipurpose SEM with smart analytical port design can accommodate not only multiple EDS but also EBSD and WDS etc. This JSM-IT510 can also be equipped with a LaB₆ electron gun, for applications requiring higher brightness and longer source lifetimes.

- Resolution: High Vacuum Mode: 3.0 nm (30 kV), 15.0 nm (1 kV)
2.0 nm (30 kV), 12.5 nm (1 kV) *LaB₆ electron gun
Low Vacuum Mode: 4.0 nm (30 kV, BED)
3.0 nm (30 kV, BED) *LaB₆ electron gun
- Landing voltage: 0.3 to 30 kV
- Direct magnification: ×5 to ×300,000 (print size of 128 × 96 mm)
- Display magnification: ×14 to ×839,724 (display size of 358 × 269 mm)



JSM-IT210

Scanning Electron Microscope

- Guided operation for necessary data acquisition -

The JSM-IT210 series is our compact thermionic-emission SEM that delivers exceptional performance at a great value. The JSM-IT210 allows for seamless navigation from optical to SEM imaging, Live EDS, Live 3D and auto functions from alignment to focus for fast, clear, and sharp images. Productivity enhancing automation is included to simplify workflow so that you can focus on results.

- Resolution: High vacuum mode: 3.0 nm (30 kV), 15.0 nm (1 kV)
Low vacuum mode: 4.0 nm (30 kV, BED)
- Landing voltage: 0.3 to 30 kV
- Direct magnification: ×5 to ×300,000 (defined with a print size of 128 × 96 mm)
- Display magnification: ×15 to ×883,652 (full screen display by standard monitor)



JCM-7000Plus NeoScope™

Benchtop Scanning Electron Microscope

The JCM-7000Plus is a benchtop scanning electron microscope with guided workflow and advanced auto functions that provides easy imaging and analysis for users at all experience levels. This versatile benchtop SEM includes high vacuum (HV) mode for observation of sub-micron structures and low vacuum (LV) for observation of non-conductive specimens in their native state. The JCM-7000 is now equipped with no code automation popular with our higher-end SEM models for enhanced efficiency and throughput.

- Resolution: High vacuum mode: 6.0 nm (20 kV)
- Direct magnification: ×10 to ×100,000 (Magnification is defined by 128 × 96 mm)
- Monitor magnification: ×24 to ×202,168 (Magnification is defined by 280 × 210 mm)
- Specimen size: 80 mm diameter
- Main options: EDS (energy dispersive X-ray spectrometer),
Stage Navigation System, Tilting and Rotating Motor Drive Holder



Ion Beam Application Equipment



JIB-PS500i FIB-SEM System

The JIB-PS500i combines high-resolution SEM imaging with a high-performance FIB to enable precise, site-specific TEM specimen preparation. A large specimen stage with wide tilt capability simplifies workflows such as block trimming, while the Check & Go function—featuring a newly developed STEM detector—allows rapid screening of prepared samples. With compatibility for double-tilt cartridges and TEM holders, the system ensures seamless correlation between FIB processing and TEM analysis.

- SEM resolution: 0.7 nm at 15 kV, 1.4 nm at 1 kV, 1.0 nm at 1 kV in BD mode
- FIB resolution: 3 nm at 30 kV
- FIB probe current: 1 pA to 100 nA



JIB-4700F MultiBeam System

The JIB-4700F features a hybrid conical objective lens, GENTLEBEAM™ (GB) mode and an in-lens detector system to deliver a guaranteed resolution of 1.6 nm at a low landing voltage of 1 kV. Using an In-Lens Schottky FEG that produces an electron beam with a maximum 300 nA probe current, this new FIB-SEM allows for high-resolution imaging and fast analyses. For the FIB column, a high-current density Ga ion beam with a maximum 90 nA probe current is employed for fast ion milling and processing of specimens.

- SEM image resolution: 1.2 nm (15 kV, GB mode), 1.6 nm (1 kV, GB mode)
- FIB image resolution: 4.0 nm (30 kV)
- FIB probe current: 1 pA to 90 nA



IB-19540CP CROSS SECTION POLISHER™ (CP)

The IB-19540CP uses a broad argon ion beam milling source to produce pristine, artifact free cross sections. Milling rates as high as 1.2 mm/h (Si equivalent, edge distance 100 μm) can be achieved with its high throughput ion source. It is possible to perform planar surface milling over a wide area up to ~20 mm in diameter with the new rotation holder option. By connecting the system to a LAN, users can check the milling status over a web browser.



IB-19550CCP COOLING CROSS SECTION POLISHER™ (CP)

The IB-19550CCP can employ a cooling function where beam sensitive specimens are cooled and maintained at a low temperature (-120 to 0 °C) while being exposed to the ion beam, thus minimizing heat damage. The unique design of the cooling function reduces the consumption of liquid nitrogen and allows for long milling times. Cooling and warming functions are automatic streamlining workflow.

In addition, a specimen can be set or removed even while the system has liquid nitrogen in the reservoir.



Peripheral Equipment



DrySD™
Energy dispersive
X-ray spectrometer



EC-52000IC
Ion Cleaner




SXES Series
Soft X-ray Emission
Spectrometer Series



JII-29080DMS
Dry Pumping Multi Station



JEC-3000FC
Auto Fine Coater

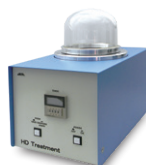

DII-29020HD
HD Treatment



EC-32010CC
Carbon Coater




DII-29010SCTR/DII-29030SCTR
Smart Coater



SM-92100EUVC
Excimer UV Cleaner



miXcroscopy™
Linked Optical and Scanning
Electron Microscope System

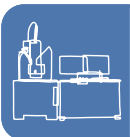

Instruments for Microarea and Surface Analysis



JXA-iHP200F

Electron Probe Microanalyzer

The JXA-iHP200F represents the next generation of Electron Probe Microanalyzer (EPMA) technology. Its Schottky field emission electron source delivers exceptional imaging and analytical resolution, combined with high and stable probe current for ultimate analytical performance. A flexible, fully customizable configuration of spectrometers and crystals allows the system to be precisely tailored to specific application requirements. The JXA-iHP200F supports panchromatic or fully quantitative hyperspectral cathodoluminescence (CL), and can be equipped with JEOL's Soft X-ray Emission Spectrometer (SXES) for advanced chemical state analysis and ultra-light element detection.



JXA-iSP100

Electron Probe Microanalyzer

The JXA-iSP100 is a highly versatile Electron Probe Microanalyzer (EPMA) that delivers high-resolution imaging and reliable analytical performance, supported by a high and stable probe current. Designed for flexibility, the system can be configured with either a tungsten (W) or LaB_6 electron source to suit a wide range of applications. A streamlined workflow, combined with customizable spectrometer and crystal configurations, enables efficient optimization for diverse analytical needs. The JXA-iSP100 supports panchromatic or fully quantitative hyperspectral cathodoluminescence (CL), and can also be equipped with JEOL's Soft X-ray Emission Spectrometer (SXES) for chemical state analysis and ultra-light element detection.



JAMP-9510F

Field Emission Auger Microprobe

The JAMP-9510F is a high-performance field-emission Auger microprobe designed for advanced surface analysis. It delivers high-throughput chemical-state characterization through a hemispherical electrostatic energy analyzer (HSA), while maintaining large probe currents even at small probe diameters enabled by its field emission gun (FEG). A eucentric tilt stage combined with an integrated charge neutralization gun allows reliable analysis of both conductive and insulating materials. With its capability to provide detailed chemical composition and chemical-state information, the JAMP-9510F offers comprehensive surface characterization for a wide range of samples.



JPS-9030

Photoelectron Spectrometer (XPS)

The JPS-9030 is a multi-purpose XPS adopting newly-developed software for greater ease-of-use. A new Kaufman-type etching ion source is installed in the specimen exchange chamber to prevent contamination of the measurement chamber. In addition to the standard Mg/Al twin anode, an infrared heating system and an Ar gas cluster ion source are available.



X-ray Fluorescence Spectrometers



JSX-1000S

Energy-dispersive X-ray fluorescence spectrometer

The JSX-1000S is an energy dispersive X-ray fluorescence spectrometer that provides high-sensitivity and high-throughput analysis across its entire energy range and can employ up to nine filters. Features include built-in recipes which automate analyses based on specific applications or sample categories (ex. ROHS), Smart FP method for accurate standard-less quantitative analysis and helium-free light element analysis in liquid by utilizing our Low Vacuum Capsule option.




Electron Diffractometer




XtaLAB Synergy-ED

Electron Diffractometer


The XtaLAB Synergy-ED is a new and fully integrated electron diffractometer that uses electron beams to create a seamless workflow from data collection to structure determination of three-dimensional molecular structures. The XtaLAB Synergy-ED is the result of combining core technologies to achieve synergy between Rigaku's high-speed, high-sensitivity detector (HyPix-ED) and single crystal analysis software (CrysAlis^{Pro} for ED) and JEOL's transmission electron microscope technology. By creating a seamless workflow from sample selection (nanocrystals) to data collection and structure analysis, the XtaLAB Synergy-ED makes electron diffraction easily accessible to non-experts who do not have the expertise in electron microscopy and crystallography that is typically required.

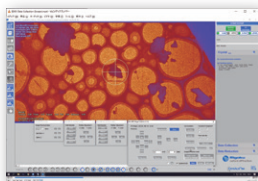

1 Select the target crystal

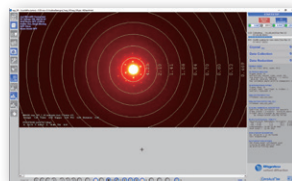


2 Measure and analyze using electron diffraction

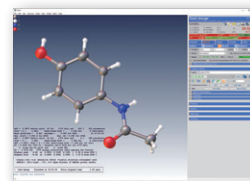

3 Analyze crystal structure



Mass Spectrometers

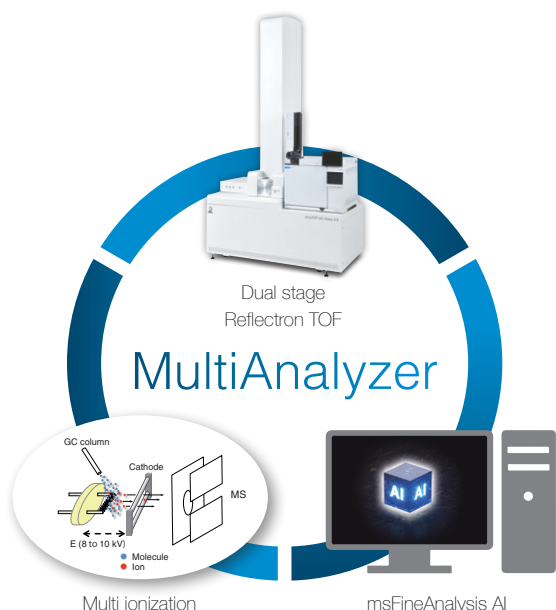

JMS-T2000GC MultiAnalyzer

Multi-ionization Unknown Compounds Analysis System


Materials are getting more complex, our environment is diversifying, and hundreds of new chemicals are created every day. The Multi-ionization Unknown Compounds Analysis System MultiAnalyzer performs structure analysis even for unknown compounds and provides answers to such complicated and challenging problems. We are committed to provide high throughput, high quality analysis results with MultiAnalyzer.



JMS-Q1600GC UltraQuad™ 2.0

Gas Chromatograph - Quadrupole Mass Spectrometer (GC-QMS)

UltraQuad™ 2.0 integrates JEOL's long-developed QMS technologies to deliver high performance, versatility, and analytical expandability across a wide range of GC-MS applications. With the standard EI ion source, it achieves an IDL of < 5 fg OFN, enabling highly sensitive measurements. Users can also select EI/PI or EI/CI combination ion sources, both of which allow ionization-method switching without venting, greatly expanding the range of analyzable compounds. The qualitative analysis software msFineAnalysis IQ uses AI-based structural estimation to identify difficult-to-analyze unknown compounds. The system also supports nitrogen carrier gas, making it a single instrument capable of meeting diverse quantitative and qualitative analytical needs.

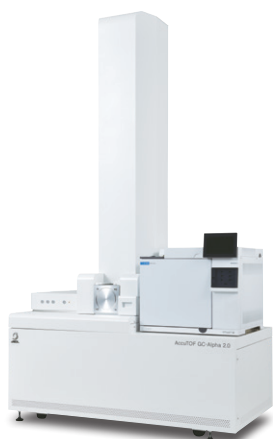

JMS-T2000GC AccuTOF™ GC-Alpha 2.0

Gas Chromatograph - Time-of-Flight Mass Spectrometer

The JMS-T2000GC AccuTOF™ GC-Alpha 2.0 is a high-performance GC-MS system that features high mass-resolution, high mass-accuracy, high sensitivity, high speed data acquisition, wide dynamic range, and wide mass range. Two types of multi-ionization combination ion sources - EI/FI/FD combination ion source and EI/PI combination ion source - are available. Soft ionization methods can produce molecular ions from unknown compounds whose molecular ions can not be produced by EI, and facilitate the elucidation of their elemental compositions. The system comes with the automatic structure analysis software "msFineAnalysis AI" for integrated structure analysis of the data measured by the multi-ionization combination ion sources, thus providing fast automatic structure analysis, which does not solely depend on mass spectral library databases.



JMS-TQ4000GC UltraQuad™ TQ

Gas chromatograph - Triple Quadrupole Mass Spectrometer (GC-QMSMS)

The new GC triple quadrupole mass spectrometer that provides both high productivity and highly sensitive analysis with the combination of JEOL's proprietary short collision cell technology and the ion ejection technology (ejecting ions after accumulating and pulsing them). It exerts its capabilities in the quantitative analysis of trace target compounds in high matrix samples, such as pesticides residues in foods. Further trace analysis is possible by using an optional ion source such as the Enhanced Performance Ion Source (EPIS).



Photo: JMS-Q1600GC + Headspace autosampler system

MS-62071STRAP

Headspace Autosampler

The MS-62071STRAP is a new-generation headspace sampler (HS) employing sorptive trap and providing ultra-high sensitivity measurement, which was not possible with the conventional HS sampler using a sample loop. In addition, the HS-GCMS system, combined with a JMS-Q1600GC, guarantees detection of mold odor in water down to 1 ppt.



JMS-S3000 NewSpiralTOF™

Ultra-High Mass-Resolution MALDI-TOFMS System

The JMS-S3000 NewSpiralTOF™ is a MALDI-TOFMS* featuring innovative SpiralTOF ion optics. Upgraded to NewSpiralTOF™, it further extends the mass range in SpiralTOF mode and enhances mass spectrometry imaging speed. Setting a new benchmark in MALDI-TOFMS performance, the JMS-S3000 delivers state-of-the-art analytical solutions for diverse research fields, including functional synthetic polymers, materials science, and biomolecules.

* matrix-assisted laser desorption/ionization time-of-flight mass spectrometer



JMS-T100LP AccuTOF™ LC-Express

Atmospheric Pressure Ionization High-Resolution Time-of-Flight Mass Spectrometer

The JMS-T100LP is a robust, easy-to-maintain, high-throughput atmospheric pressure ionization time-of-flight mass spectrometer intended for high productivity with multiple ionization methods. The AccuTOF™ LC-Express is the advanced version of JMS-T100LP with a significantly improved dynamic range. By using Direct Analysis in Real Time (DART™), which is JEOL's proprietary ionization technology of rapidly providing accurate mass information, as well as easy replacement of the ion source with the electrospray ionization (ESI) source for LC-MS operation or with the ColdSpray ionization (CSI) source, the AccuTOF™ LC-Express satisfies various research needs in organic chemistry, material science, and other fields.



JMS-T100LP + DART™ Ion Source

Ambient Ionization High-Resolution TOFMS

Direct Analysis in Real Time (DART™) is an ion source that can analyze samples with various shapes and states without any sample preparation. DART™ was born in 2003 at the mass spectrometry applications laboratory of JEOL USA, Inc. Among a series of new ionization techniques, which were later termed "ambient ionization," DART™ was the first to have been invented and the first to have been commercialized in 2005. It is possible to acquire high mass-resolution, accurate-mass spectra in real time by simply presenting samples of various shapes and states to the DART™ ion source without any sample preparation. DART™ can handle samples with arbitrary shapes or "dirty" samples that conventional analytical methods struggle. DART™ was developed for the JEOL AccuTOF™ series of mass spectrometers. AccuTOF™ LC-Express and DART™ are the perfect combination.



NMR Spectrometers



JNM-ECZL (ECZ Luminous™) series FT NMR System

The JNM-ECZL series (ECZ Luminous™) is an FT NMR system equipped with state-of-the-art digital and high-frequency technologies. The highly integrated Smart Transceiver System, a high-speed, high-precision digital high-frequency control circuit, enables further miniaturization and ensures the stability of spectrometer. The new Multi Frequency Drive System enables multi-resonance measurements in a standard configuration, providing a wider range of solutions.

ECZL G series

This is a flagship model that supports diversified, cutting-edge applications, including high-magnetic-field, and solid-state NMR measurements. The system is flexible in terms of expansion, with support for three or more channels, high-power amplifiers, and high-power magnetic field gradients.

- Frequency: 400 MHz to 1.3 GHz
- Sample type: solution/solid
- Number of channels: 2 (standard), expandable up to 8 channels
- Power amplifier: High frequency side
400 to 600 MHz, 100 W, 700 MHz or higher, 200 W
500 W/1000 W (as an option)
Low frequency side
400 to 600 MHz, 300 W, 700 MHz or higher, 500 W
1000 W/2000 W (as an option)
- Magnetic field gradient amplifier: 10 A (standard), 30 A/50 A (as an option)

ECZL R series

This is a compact model that supports solid-state NMR measurements. Its footprint is less than 50% of the spectrometer of JNM-ECZR.

- Frequency: 400 MHz to 600 MHz
- Sample type: solution/solid
- Number of channels: 2
- Power amplifier: High frequency side 100 W, Low frequency side 300 W
- Magnetic field gradient amplifier: 10 A

ECZL S series

This is an entry-level model dedicated to 400 MHz solution NMR while incorporating the high-performance digital high-frequency technology of the ECZL series.

- Frequency: 400 MHz
- Sample type: solution
- Number of channels: 2
- Power amplifier: High frequency side 50 W, Low frequency side 150 W
- Magnetic field gradient amplifier: 10 A



Cryogen Reclamation System (CR-80)



Nitrogen Replenishment System (NR-50)



Liquid Nitrogen Generator (NS-20W-J)

CR-80 Cryogen Reclamation System

The CR-80 is an attachment that can greatly reduce the cryogen evaporation of both liquid helium and liquid nitrogen which are inevitable for superconducting magnets for NMR. This system will greatly help in the management of cryogens for NMR system by eliminating concerns about follow-up services as it is an in-house product, concerns about unstable supply of cryogens, price increases, and scheduling of filling operations. We also have the "NR-50" which is a system to reduce evaporation of liquid nitrogen and the "NS-20W-J" which is an auto supplier of liquid nitrogen allowing for liquid nitrogen to be taken out.



Solution NMR probe

Solution NMR probe

In solution NMR, samples dissolved in a solvent are measured by placing them in glass tubes. A wide variety of probes can be provided depending on the magnet size, target nuclei, and sample tube size. The newly developed cooling type probe, SuperCOOL MARVEL probe offers a dramatic improvement in sensitivity compared with conventional probes of the same class. It is user-friendly broadband type, enabling high sensitivity measurements across multiple nuclei. This probe is suitable not only for routine measurements, but also for advanced application experiments.



SuperCOOL MARVEL probe




Solid-state NMR probe

Solid-state NMR accessories

In solid-state NMR, solid-state samples such as powders and films are placed in a dedicated ceramics sample tube for measurement. The solid-state NMR technique is effective for analyzing samples that do not dissolve in solvents such as inorganic materials, and crystalline polymorph that can be meaningfully measured in the solid state. Use of AUTOMAS probe, ROTORCARRIER™, and auto sample changer enables automatic measurement of multiple samples like solution NMR.



ROTORCARRIER™ for AUTOMAS probe

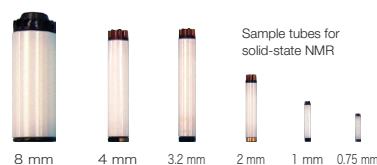

Sample tubes for solid-state NMR



ESR Spectrometers




JES-X3 series ESR Spectrometers

Recently, it has been widely accepted that even relatively few unpaired electrons in a sample can affect the function of the material, so a lower detection limit (higher sensitivity) is required of ESR measurements. The JES-X3 Series has achieved higher sensitivity by developing a low-noise Gunn oscillator for its new spectrometer.

	X310	X320
Maximum Magnetic Field	0.65 T	1.3 T
Sweep Width	±0.01 to 250 mT	±0.01 to 500 mT
Pole Gap	60 mm	60 mm
Frequency Range	8.750 to 9.650 GHz	
Field Resolution	2.35 μ T	
Correction by Marker	Standard	



Semiconductor Equipment

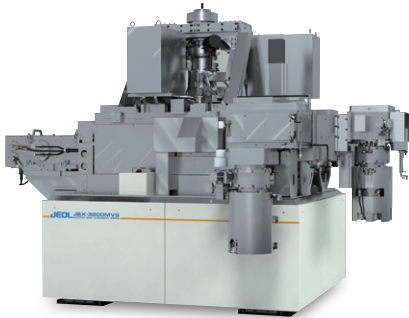


JBX-3200MVS

Electron Beam Lithography System

The JBX-3200MVS is a variable-shaped electron beam lithography system for mask making of 32 nm to 28 nm nodes. Its advanced technology achieves high speed, high precision and high reliability. This EB system uses a variable-shaped 50 kV electron beam and a step-and-repeat specimen stage.

- High throughput enabled by a new PEC (Proximity Effect Correction) system
- High-speed data-transfer system for a large quantity data
- High-precision design for mask/reticle production

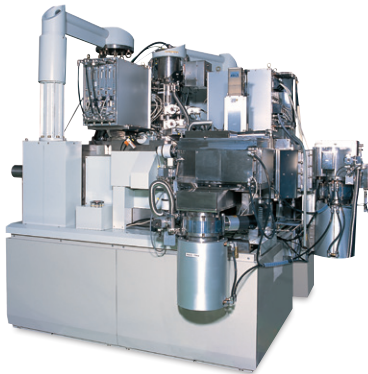

JBX-3050MV/S

Electron Beam Lithography System

The JBX-3050MV/S is a variable-shaped electron beam lithography system for mask making up to the 45 nm nodes. Its advanced technology achieves high speed, high precision and high reliability. This EB system uses a variable-shaped 50 kV electron beam and a step-and-repeat stage.



JBX-A9

Electron Beam Lithography System

The JBX-A9 is a spot-beam e-beam lithography system for 300 mm wafers that inherits the basic specifications of the JBX-9500FS while achieving power and space savings and is environmentally friendly through the use of a refrigerant-free chiller. This system is suitable for applications that require particularly high beam positioning accuracy, such as photonic crystal device fabrication.

- Accelerating voltage: 100 kV
- Electron beam source: ZrO/W Schottky emitter
- Scanning speed: 200 MHz to 250 Hz
- Substrate size: Maximum 300 mm wafer



JBX-8100FS

Electron Beam Lithography System

The JBX-8100FS is a spot-beam e-beam lithography system that performs superior accuracy and writing speed. Upgradable platform offers various options to meet the fields from high-end nanostructure applications to batch production of compound semiconductor.

- Accelerating voltage: 100 kV (option: 200 kV/130 kV/50 kV/25 kV)
- Electron beam source: ZrO/W Schottky emitter
- Scanning speed: 125 MHz to 250 Hz
- Substrate size: Maximum 200 mm wafer



Metal AM Machine



JAM-5200EBM


Electron Beam Metal 3D Printer

Metal 3D printer that employs additive manufacturing with powder bed fusion techniques. This 3D printer irradiates the laid metal powders with an electron beam to melt and solidify the powders for fine-structure production. The features of the printer include direct manufacturing of 3D micro-structures, and large-output, high-speed, high-density manufacturing, thus allowing for additive manufacturing with high quality and high reproducibility.


- Manufacturing dimensions: Maximum 250 mm (dia.) × 400 mm (H)
*Depends on the building material and specifications
- Electron beam output: Maximum 6 kW




Industrial Equipment for thin-film formation and material processing

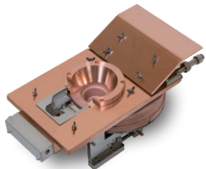


BS-60211DEM

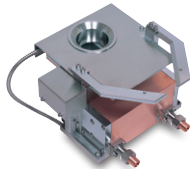

BS-60060DEBS

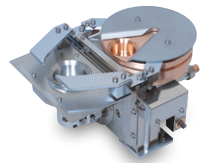

EBG-102UB6S

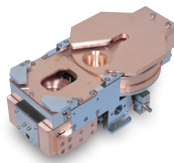

EBG-203UB6S


BS series/EBG series

Electron Beam Source

Electron beam sources for vacuum evaporation of thin-film deposition. Various models are selectable for oxide films, nitride films and metal films. Extended-use models are also available for multilayer films and continuous evaporation with the long-lifetime filament.



BS-60610BDS

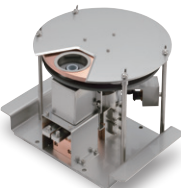

BS-60310BDS

BS-60610BDS/BS-60310BDS

Bombardment Deposition Source

Deposition sources that utilize electron-beam bombardment indirect heating techniques. Suitable for low-damage, low-defect, low-absorption film evaporation. The BS-60610BDS achieves large capacity and high rate for liners.



Industrial Equipment for thin-film formation and material processing



JEBG-3000UB

JEBG-1000UB

JEBG series

High-power Electron Beam Source

High-power electron beam sources for vacuum evaporation of metal and metal-oxide for wide plastic films or large glass plates that are continuously fed. They can also be used for vacuum melting of high-melting point metals.



BS-80011BPG

BS-80020CPPS

BS-80011BPG/BS-80020CPPS

Plasma Source

Plasma Sources are installed in a vacuum chamber and generate high-density plasma. Used for Ion Plating (Plasma Assisted Deposition) and it is possible to improve film properties for optical thin films, protective films and functional films. Because high density plasma can be generated in a mass space, high-rate deposition to a large area is possible.



BS-04 series

Rotary Sensor

Multi-point sensor for measuring the deposition rate and film thickness in vacuum evaporation or sputtering by connecting it to a quartz-rate control monitor. Equipped with 6 or 12 quartz crystals, smooth crystal replacement is achieved for deposition of multilayer or thick films. Sensor length in the chamber is selectable and two kinds of sensor-head angle are available, thus enabling flexible sensor installation into small to large equipment.



BS series

Electron Beam Source for research and development

Electron beam source suitable for research and development. Various optional accessories can be added for a wide range of customization, including lift-off evaporation with low temperature, low damage, low defect, etc.



TP-40020NPS Nano-powder synthesis system

TP series

RF Induction Thermal Plasma System

Thermal plasma around 10,000 °C by RF inductive coupling. Applications of the thermal plasma system include nano-powder synthesis, fine-powder spheroidization, high-speed thick-film synthesis and gas decomposition. Irrespective of materials and gases, chemical reactions (oxidization, nitriding) and surface reforming can be made.



TP-99260FDR

Powder Feeder

The table-type powder feeder that can stably feed minute particles of various sizes ranging from submicron to 100 μm or larger. It feeds particles using carrier gas through the piping. Even particles of a few microns or less in size that are highly cohesive and therefore poor in fluidity can be fed continuously. The powder feeding rate is approximately from 0.1 to 100 g/min. The closed structure enables feeding into a reduced-pressure atmosphere. Not only single-composition particles but also mixed-composition particles can be fed continuously without changing the mixing ratio. It also supports real-time monitoring of powder feeding, feedback, and external control.



Certain products in this brochure are controlled under the "Foreign Exchange and Foreign Trade Law" of Japan in compliance with international security export control. JEOL Ltd. must provide the Japanese Government with "End-user's Statement of Assurance" and "End-use Certificate" in order to obtain the export license needed for export from Japan. If the product to be exported is in this category, the end user will be asked to fill in these certificate forms.



3-1-2 Musashino Akishima Tokyo 196-8558 Japan Sales Division Tel. 81-3-6262-3561 Fax. 81-3-6262-3577
www.jeol.com ISO 9001 • ISO 14001 Certified

ARGENTINA
COASIN S.A.C.I.F.
Virrey del Pino 4071,
C1430CAM-Buenos Aires
Argentina
Tel. 54-11-4552-3185
Fax. 54-11-4555-3321

AUSTRALIA & NEW ZEALAND
JEOL (AUSTRALASIA) Pty.Ltd.
Suite 1, Level 2, 18 Aquatic Drive
Frenchs Forest NSW 2086
Australia
Tel. 61-2-9451-3855
Fax. 61-2-9451-3822

AUSTRIA
JEOL (GERMANY) GmbH
Gute Aenger 30
85356 Freising, Germany
Tel. 49-8161-9845-0
Fax. 49-8161-9845-100

BANGLADESH
A.Q. CHOWDHURY SCIENCE & SYNERGY PVT. LTD.
87, Suhrawardy Avenue, Floor 2
Baridhara, Dhaka 1212
Bangladesh
Tel. 880-2-2-2226-2272
Fax. 880-2-2-2226-4428

BELGIUM
JEOL (EUROPE) B.V.
Planet II, Gebouw B
Leuvensesteenweg 542,
B-1930 Zaventem
Belgium
Tel. 32-2-720-0560
Fax. 32-2-720-6134

BRAZIL
JEOL Brasil Instrumentos Cientificos Ltda.
Av. Jabaquara, 2958 5° andar conjunto 52
04046-500 Sao Paulo, SP
Brazil
Tel. 55-11-5070-4000
Fax. 55-11-5070-4010

CANADA
JEOL CANADA, INC.
3275 1ere Rue, Local #8
St-Hubert, QC J3Y-8Y6, Canada
Tel. 1-450-765-1885
Fax. 1-450-766-1987

CHINA
JEOL (BEIJING) CO., LTD.
Zhongkeziyuan Building South Tower 2F,
Zhongguancun Nanshanjie Street No. 6,
Haidian District, Beijing 100190, P.R.China
Tel. 86-10-6804-6332/6322/6323
Fax. 86-10-6804-6324

JEOL (BEIJING) CO., LTD., SHANGHAI BRANCH
2F-BC Room, Building A, Mingji Business Plaza,
No.207 Songhong Road, Changning District,
Shanghai 200335, P.R.China
Tel. 86-21-6249-4487
Tel. 86-21-5836-8350
Fax. 86-21-5836-3668

JEOL (BEIJING) CO., LTD., GUANGZHOU BRANCH
Rm.3501, OneLinkCenter, 230 Tianhe Road, Tianhe District,
Guangzhou, Guangdong Prov., 510620, China
Tel. 86-20-8778-7848
Fax. 86-21-8778-4268

JEOL (BEIJING) CO., LTD., WUHAN BRANCH
R1104, Chicony Center, No.10 Luoyu Road,
Hongshan District, Wuhan, 430070, P.R.China
Tel. 86-27-8713-2567
Fax. 86-27-8713-2567

JEOL LTD. (BEIJING) CO., LTD., CHENGDU BRANCH
1807A Zongfu Building,
NO. 35 Zhongfu Road, Chengdu, Sichuan, 610016
P.R. China
Tel. 86-28-8662-2554
Fax. 86-28-8662-2564

EGYPT
JEOL SERVICE BUREAU
3rd Fl. Nile Center Bldg., Nawal Street,
Dokki, (Cairo), Egypt
Tel. 20-2-3335-7220
Fax. 20-2-3338-4186

FRANCE
JEOL (EUROPE) SAS
Espace Claude Monet, 1 Allée de Giverny
78290, Croissy-sur-Seine, France
Tel. 33-13015-3737
Fax. 33-13015-3717

GERMANY
JEOL (GERMANY) GmbH
Gute Aenger 30
85356 Freising, Germany
Tel. 49-8161-9845-0
Fax. 49-8161-9845-100

UNITED KINGDOM (UK)
JEOL(UK)LTD.
Silver Court, Watchmead,
Welwyn Garden City, Hertfordshire AL7 1LT, UK
Tel. 44-1707-37-7117
Fax. 44-1707-37-3254

GREECE
N. ASTERJADIS S.A.
56-58 S. Tirkoup Stn, P.O. Box 26140
GR-10022, Athens, Greece
Tel. 30-1-823-5383
Fax. 30-1-823-9567

HONG KONG
FARMING LTD.
Unit No. 1009, 10/F, Prosperity
Millennia Plaza, 663 King's Road, North Point, Hong Kong
Tel. 852-2815-7299
Fax. 852-2581-4635

TechMax (Hong Kong) Company Limited
Sales & Service (NMR/ESR/MS)
16/F, 700 Nathan Road,
Mongkok, Kowloon, Hong Kong
Tel. 852-5890-7016
Fax. 852-3102-2805

INDIA
JEOL INDIA PVT. LTD.
Unit No.305, 3rd Floor,
ABW Elegance Tower,
Jasola District Centre,
New Delhi 110 025, India
Tel. 91-11-4595-8000/8005/8017

JEOL INDIA PVT. LTD. Mumbai Office
Business Suite No. S31,S32, S33
Vatika Business Centre - Unit No. G5, Ground Floor and
Unit No. 02, First Floor
Trade Centre, Bandra Kurla Complex
Bandra (East), Mumbai - 400051, INDIA
Tel. 91-76-6697-7120

JEOL INDIA PVT. LTD. JEOL SOUTHY CENTRE (JSOL)
Anthill IQ Workspaces, Unit No. G-01, No. 650/A,
Kodigehalli Extension, Bellary Road, Hebbal,
Bengaluru, 560 092, India
Tel. 91-87964-42640

JEOL INDIA PVT. LTD. Chennai Office
Door No.923, 9th Floor, Regus Olympia Platina
Plot No.33-B, South Phase, Guindy Industrial Estate,
Chennai - 600032, India
Tel. 91-99-4009-7547

JEOL INDIA PVT. LTD. Kolkata Office
Office Number (CM DAL'HOUSE)
Regus, PS Arcadia, 9th Floor, 4A, Camac Street,
Abanindra Nath Thakur
Sarani, Kolkata, West Bengal- 700016, India
Tel. 91-98-3023-0484

JEOL INDIA PVT. LTD. Hyderabad Office
433, Regus Solitaire Business Centre,
1-10-39 to 44, Level 4,
Gundlupet Towers, Old Airport Road,
Begumpet, Hyderabad - 500016, India
Tel. 91-97-4111-5885

INDONESIA
PT. TEKNO LABINDO Penta Perkasa
Komplek Gading Bukit Indah Blok I/11
Jl. Bukit Gading Raya Kelapa Gading Permai,
Jakarta 14240, Indonesia
Tel. 62-21-4584-7057/7058
Fax. 62-21-4584-7279

ITALY
JEOL (ITALIA) S.p.A.
Palazzo Piacinotti - Milano 3 City,
Via Ludovico il Moro, 6/A
20079 Bassiglio(MI), Italy
Tel. 39-02-904-1431
Fax. 39-02-9041-4343

KOREA
JEOL KOREA LTD.
AJ Bldg. 6F,
9, Jeongui-ro 8-gil, Songpa-gu, Seoul,
05836, Korea
Tel. 82-2-511-5501
Fax. 82-2-511-2635

KUWAIT
Ashraf & CO. Ltd.,
Block 2 Bldg #277, Ardya Industrial Area,
P.O. Box 3555 13036 Safat, Kuwait
Tel. 965-2209-0750

MALAYSIA
JEOL (MALAYSIA) SDN.BHD.
508, Block A, Level 5,
Kelana Business Center,
97, Jalan SS 7/2, Kelana Jaya,
47301 Petaling Jaya, Selangor, Malaysia
Tel. 60-3-7492-7722
Fax. 60-3-7492-7723

MEXICO
JEOL DE MEXICO S.A. DE C.V.
Arkansas 11 Piso 2
Colonia Napoléon
Delegación Benito Juárez, C.P.03810
Mexico D.F., Mexico
Tel. 52-55-5448-5900
Fax. 52-55-5211-0720

UNITED ARAB EMIRATES (UAE)
JEOL GULF FZCO
P.O. Box No. 371107
Dubai Airport Free Trade Zone West Wing SWA No. G12,
Dubai, UAE
Tel. 971-4-609-1497
Fax. 971-4-609-1498

PAKISTAN
ANALYTICAL MEASURING SYSTEM (PVT) LTD. (AMS LTD.)
14-C Main Sehar Commercial Avenue Lane 4,
Khayaban-e-Sehar,
D.H.A-VII, Karachi-75500, Pakistan
Tel. 92-21-3534-5581/0747
Fax. 92-21-3534-5582

PANAMA
PROMED S.A.

Parque Industrial Costa del Este
Urbanización Costa del Este
Apartado 0816-01755, Panama
Tel. 507-303-3100
Fax. 507-303-3115

PHILIPPINES
JATEC Philippines Corporation
28 Floor, The Enterprise Center Tower 2,
Ayala Avenue corner Paseo de Roxas,
Ermy, San Lorenzo, Makati City, 1226, Philippines
Tel. 63-2-8849-3904

PORTUGAL
Izasa Portugal Lda.
Av. do Forte N6 Piso 3, Sala 2,24,
2790-072 Camarnide, Portugal
Tel. 351-21-424-73-00
Fax. 351-21-418-60-20

QATAR
Mannai Trading Company W.L.L.
Al J Emadi Complex,
Salwa Road P.O.Box 76, Doha, Qatar
Tel. 974-4455-8216
Fax. 974-4455-8214

RUSSIA
JEOL (RUS) LLC
Office 351, floor 3, 23,
Novolobodskaya St,
Moscow 127055, Russia
Tel. 7-495-748-7791/7792
Fax. 7-495-748-7793

SAUDI ARABIA
ABDULREHMAN ALGOSABI G.T.C. (Riyadh)
Algosabi Building-Old Airport Road
P.O. Box 215, Riyadh-11411, Saudi Arabia
Tel. 966-11-272-6999
Fax. 966-11-477-1374

SWEDEN
JEOL (Nordic) AB
Hammarbacken 6A, SE-19127 Sollentuna
Sweden
Tel. 46-8-28-2800
Fax. 46-8-29-1647

SINGAPORE
JEOL ASIA PTE.LTD.
2 Corporation Road
#01-12 Corporation Place
619494, Singapore
Tel. 65-6565-6869
Fax. 65-6565-7552

SOUTH AFRICA
ANGSTROM Scientific (Pty) Ltd.
370 Angus Crescent, Northlands Business Park
29 Newmarket Road
Northridge, Randburg 2169, South Africa
Tel. 27-11-462-1347

SPAIN
PALEX MEDICAL S.A.
Calle Jesus Serra Santamans, Num. 5
08174 Sant Cugat Del Valles, Barcelona
Spain
Tel. 34-900-878-233

SWITZERLAND
JEOL (GERMANY) GmbH
Toldstrasse 60
8002 Zurich
Schweiz
Tel. 41-4428-36124
Fax. 41-4429-36123

TAIWAN
JIE DONG CO., LTD.
7F, 112, Chung Hsiao East Road,
Section 1, Taipei 10023, Taiwan
Tel. 886-2-2395-2978
Fax. 886-2-2392-4655

For NMR & Mass Spectrometer Products:
TechMax Technical Co., Ltd.
17F, No. 216, Sec. 3, New Taipei Blvd., Xinzhuang Dist.,
New Taipei City 242032, Taiwan
Tel. 886-2-8990-1779
Fax. 886-2-8990-2559

For Semiconductor Products:
JEOL TAIWAN SEMICONDUCTORS LTD.
No.198, Xingong 2nd Rd., Zhubei City,
Hsinchu County 302004, Taiwan
Tel. 886-3-571-5556
Fax. 886-3-571-5151

THAILAND
BEETHAI BANGKOK EQUIPMENT and CHEMICAL CO., Ltd.
99/9 Moo 2, Maha Sawat, Putthamonthon, Nakhon Pathom, 73170
Thailand
Tel. 66-3424-5299
Fax. 66-3424-5250

JEOL ASIA(THAILAND)CO., LTD.
118 BLOCK 28X Building, 4th Floor, Room No.X401, Soi Chula 5,
Wangmai sub-district, Phlathumwan district, Bangkok 10330
Thailand
Tel. 66-2121-7424

THE NETHERLANDS
JEOL (EUROPE) B.V.
Lireweg 4, NL-2153 PH Nieuw-Vennep,
The Netherlands
Tel. 31-232-623-600
Fax. 31-232-623-501

TURKEY
KARFO Endustriyel
Yukan Dudullu, Karlo Plaza, Şerifali Mah, Beyan Sokağı No:23,
34775 İstanbul, Turkey
Tel. 90-216-660-0315

META Analitik
Atatürk Mah, Meriç Cad, Turkuaz Plaza No:5/38 Ataşehir,
34758 İstanbul, Turkey
Tel. 90-216-680-6020

UNITED STATES OF AMERICA (USA)
JEOL USA, INC.
11 Dearborn Road, Peabody, MA 01960, USA
Tel. 1-978-636-5800
Fax. 1-978-536-2205/2206

JEOL USA, INC. WEST OFFICE
5653 Stoneridge Drive Suite #110
Pleasanton, CA 94588, USA
Tel. 1-925-737-1740
Fax. 1-925-737-1749

VENEZUELA
GOMSA Service and Supply C.A.
Urbanización Montalban III
Residencias Don Andres - Piso 7 - Apartamento 74
Avenida 3, entre calles 7 y 6
Montalban, Caracas, Venezuela
Tel. 58-212-443-4342
Fax. 58-212-443-4342

VIETNAM
T&N Trading and Investment Company Limited
19 Hang Thiec, Hoan Kiem Ward,
Hanoi, Vietnam
Tel. 84-94343-1990/Ho Chi Minh (HCM) 84-90565-3866